


Email :
sales@schooleducationalinstruments.com
Phone: +91-0171-2601773

Product Name :
Osborne Reynolds Demonstration

Product Code :
SCHOOL-FUD810025

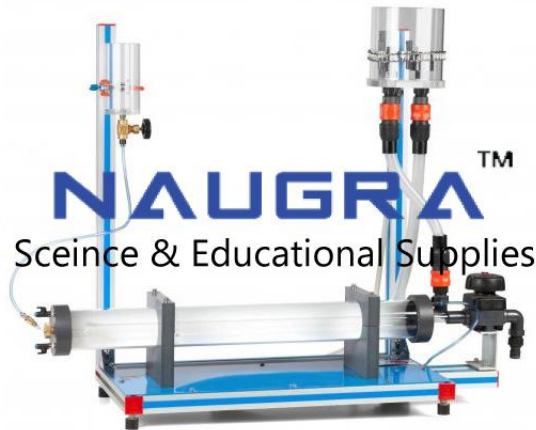

Description :

Osborne Reynolds Demonstration

Technical Specification :

The Osborne-Reynolds experiment is used to study the characteristics of a liquid flow through a pipe.

It is also used to determine the Reynolds number at each state of the liquid.

This module makes it possible to study the characteristics of the flow of a liquid inside a pipe and the behavior of such flow.

Besides, it is possible to determine the range of the laminar and turbulent flows using the Reynolds number.

Thus, difference between laminar, turbulent and transition flows can be demonstrated and the Reynolds number can be calculated for each regime.

This module consists on a transparent and horizontal pipe section, which makes it possible to visualize the fluid, a water supply tank, which guarantees the flow homogeneity, and a needle connected to a tank through a hose, from where the dye is supplied.

Water flow in the test section can be regulated by means of a valve.

Water can be supplied either using the Hydraulics Bench or from the Basic Hydraulic Feed System.

Observation of laminar, transition and turbulent flows.

Association of laminar, transition and turbulent flows with their corresponding Reynolds number.

Observation of the parabolic velocity profile.

This module is mounted on an anodized aluminum structure with painted steel panel.

Methacrylate test pipe with an airfoil-shaped inlet section:

NAUGRA[®]

School Educational Instruments

Website: [www.https://www.schooleducationalinstruments.com/](https://www.schooleducationalinstruments.com/), **Email:** sales@schooleducationalinstruments.com

Address: Ambala Cantt, Haryana, India **Phone:** 91-0171-2643080